

# Gymna Stools

---

*Take your seat*



**gymna**<sup>®</sup>

[gymna.com](http://gymna.com)



## Gymna Saddle Stool and Gymna Stool

*Enhancing the professional  
image of your practice*

Gymna Stools and Saddle Stools perfectly match with Gymna treatment tables. Together they enhance the professional image of your practice.

Saddle Stool:	Width 42 cm, with adjustable pelvic tilt
Stool:	Diameter 35 cm, with 360° lever for easy adjustment
Upholstery:	Classique (PVC) 8 colours Premium (PU) 6 colours
Height:	Gymna Stool 54-73 cm Gymna Saddle Stool 56-75 cm
Max. load:	150 kg
Base:	Aluminium with white powder coating
Wheels:	Self-locking casters, soft tread for harder surfaces

### Gymna Stool

The round and flat seat with a 35 cm diameter gives the stool a universal character. The Gymna Stool is suited for use by patients as well. Height adjustment is easy due to the 360° lever that is always accessible from all sides.



### Upholstery

The seats are elegantly finished with quality details in durable, phthalate-free and biocompatible synthetic leather: Classique or Premium. Gymna Stools are available in the same colours as the gymna.PRO series.

### Gymna Saddle Stool


This is the ideal ergonomic seat for daily use by the medical professional. Due to the saddle shape, the user sits in a more open hip angle so that the pelvis tilts forward and the spine will assume the natural S-shape. The pelvic tilt can also be adjusted separately.



- Active and improved sitting position, lowering back strain.
- Improved mobility while sitting

### Colours

#### Premium

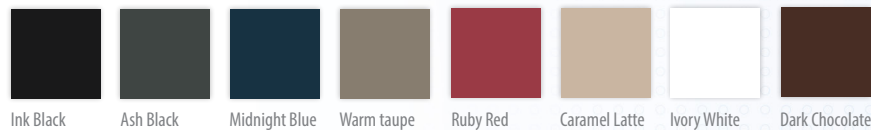

#### Classique



### Height adjustment

The stools are adjustable in height by means of a reliable gas spring. The stools are best suited for people between 160 - 185 cm. The maximum load is 150 kg.



### 5-star base

The aluminum 5-star base is coated with a scratch-resistant white powder and has a diameter of 50 cm. this is sufficient to guarantee stability yet small enough to be used effectively near the treatment table.



### Wheels

The five self-locking casters contribute to safety in practice. The stool will not suddenly slip when you sit down. The stool rolls easily as soon as you take a seat. The wheels have a soft tread, making them especially suitable for hard surfaces, commonly used in most practices.



Cryotherapy



Tecar - Diathermy



Shockwave Therapy



Physio Care



Treatment Tables




Electrolysis



Electrotherapy



Your dealer



[gymna.com](http://gymna.com)

**gymna**<sup>®</sup>

GymnaUniphy NV Pasweg 6A | B-3740 Bilzen, Belgium

Tel: +32 (0)89 510 510