

# *gymna*.ONE

*The art of therapy*

DATASHEET



***gymna***<sup>®</sup>

# *gymna.ONE – unique in its style*

*gymna.ONE is not just a treatment table—it's a reflection of unique style and professional aspirations.*

gymna.ONE offers unmatched durability, ergonomic design, and advanced technology, ensuring superior comfort and efficiency for healthcare professionals.

## **Fast and secure height adjustment**

Raising the table from lowest to highest position is significantly quicker compared to previous models, saving precious time for both therapist and patient.

**Gymna i-Control**, the safest height adjustment system, is positioned in a way that it is easily accessible from anywhere around the gymna.ONE table without the risk of unintended movement.





## Ergonomic design ensures optimal patient comfort

Designed to ensure optimal patient comfort and ergonomics gymna.ONE comes with real Gymna signature features such as our **One-Click arm support system**, **Ergomax face support** and **cushion heating**.




Lower both arm supports at the same time with the One-Click arm support system. The arm supports are slightly tilted for a more relaxed cervical and shoulder. The soft and seamless Ergomax face support is made of waterproof and washable material for optimal hygiene. Pair these two with our cushion heating and create a comfortable environment for your patient.



## Exquisite design combined with superior craftsmanship

The organically shaped aluminium arms, the exquisitely designed multi-layered wooden base and the functional but elegant edge protection, make gymna.ONE a true work of art. Combining craftsmanship with design, results in more than just a functional treatment table. It transforms any space into an inviting sanctuary of comfort and style.

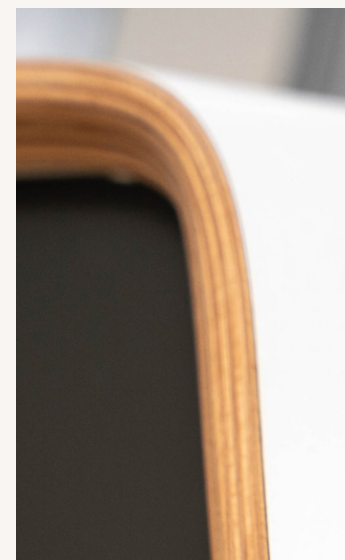


**SPECIFICATIONS**

Dimensions:	206 x 67 cm
Weight:	112 - 127 kg
Height adjustment:	50 – 100 cm
Max. lifting capacity:	200 kg
Motor force:	10.000 N
Time to lift:	18 sec.
Mains voltage:	100-240VAC, 50-60 Hz
Maximum consumption:	3,0 A
Insulation class:	II 
Device Classification:	I 
Applied parts:	Type BF
Humidity resistance:	IP21
Conformity:	EC directive 2017/745, 2012/19, 2011/65 IEC/EN 60601-1 IEC/EN 60601-1-2
Warranty:	Table 5 years, Premium upholstery 3 years

**FUNCTIONALITY**

gymna.ONE series	D6	T7	Q8
2-section table	●		
3-section table		●	
4-section table			●
Head section	●	●	●
Ergomax face support	●	●	●
One-Click arm supports	●	●	●
Fold away lateral supports	●	●	●
Cushion heating	●	●	●
Premium PU upholstery	●	●	●
Lumbar flexion, electric		●	
Back section		●	
Back section, electric		○	●
Sitting section, electric			●
Lower leg section			●
Electric height adjustment	●	●	●
i-Control switch	●	●	●
Safety switch	●	●	●
Wheel system	●	●	●
4 colour combinations	●	●	●
LED pack	○	○	○

● = standard ○ = optional

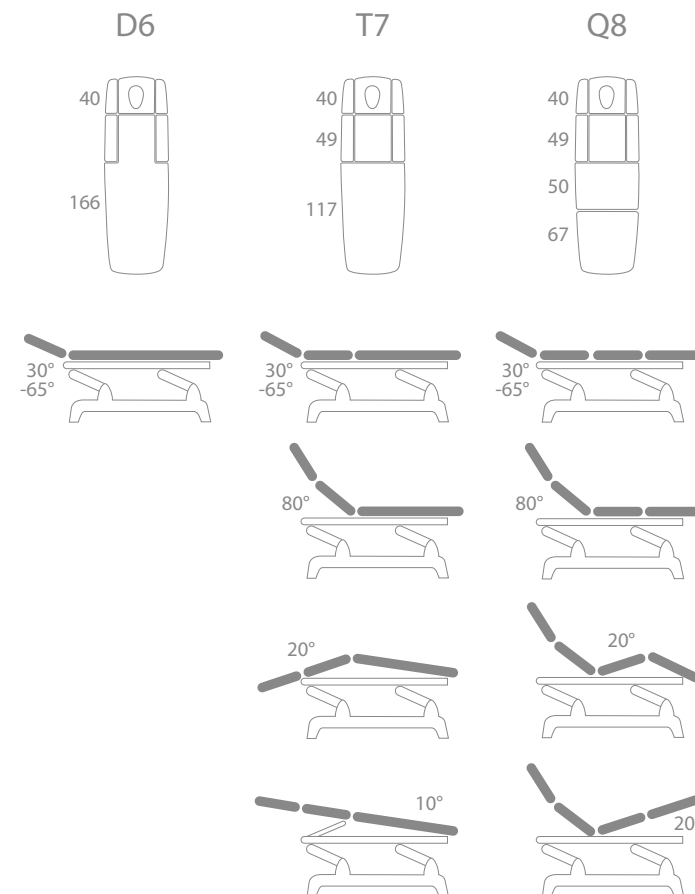**ADJUSTABILITY****CUSHIONS**

Upholstery:	PU vegan leather, phthalate free, biocompatible, abrasion resistant and flame retardant.
Foam:	PE, dual padded. Lower layer 100 kg/m <sup>3</sup> Upper layer 40 kg/m <sup>3</sup>
Cushion heating:	34°C, 36°C and 38°C Bluetooth wireless operation Automatic timer 30 minutes Max. consumption = 180W
Ergomax:	PU foam, closed skin, low density (5° Shore)

**CUSHION COLOURS**

Midnight Blue	Ink Black	Ash Black
Ruby Red	Ivory White	Caramel Latte
Warm taupe	Dark Chocolate	Cognac Brown

gymna.com ***gymna***<sup>®</sup>

GymnaUniphy NV Pasweg 6A | B-3740 Bilzen, Belgium

Tel: +32 (0)89 510 532