

Product datasheet

gymna[®]

Pulson 400M

The Pulson 400M, for ultrasound. Including a dedicated Mobile 400.

The Pulson 400M, touchscreen device for ultrasound applications. Equipped with multifrequency ultrasound heads for superficial and deeper lying pathologies. Both ultrasound heads with 1 cm² and 4 cm² are standard accessories.

Standard supplied with a functional Mobile 400 allowing easy storage and access to all accessories.

The color touchscreen, for intuitive navigation, guides you to optimal treatment protocols and access to extended clinical information with GTS.



GUIDED THERAPY SYSTEM



Article number: 360 452 Pulson 400M white
360 652 Pulson 400M black

Characteristics

■ Ultrasound:

- Multifrequency heads (1 and 3 MHz), 1 cm² and 4 cm²
- Continuous and pulsed mode (10–20–30–40–50–100%)
- Acoustic and visual contact control led
- 2 ultrasound output connectors

■ Functionalities:

- Objectives: 15
- Indications: 23
- Diagnostic: 1
- Body Area: 23
- Default therapy programs: 2
- Free memory: 550
- Anatomical library: 91

Technical specifications

Languages : 13

Mains voltages : 100-240-VAC, 50/60 Hz +/- 10%

Max. Power—in operation : 100 VA

Device (b x h x d) : 360 x 260 x 285 mm

Weight incl. accessories : c.a 7,8 kg

Mobile 400 (b x h x d) : 500 x 1100 x 600 mm

Weight Mobile : c.a 36 kg

Electrical safety protection : Class II

Applied parts : Type BF

MDD Classification : IIa

Conformity : Directive MDD 93/42/ECC

Standard accessories

100 689	Mains lead
360 111	US-head 1&3 MHz 1cm ² US 401
360 114	US-head 1&3 MHz 4cm ² US 404
341 088	Contact gel, 500 ml
340 505	Touch pen
360 808	Gymna Mobile 400



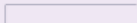



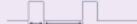







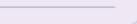



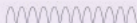

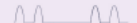

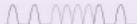

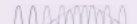












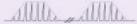
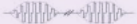
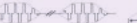




Manuals



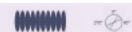


323 011	Safety Instructions
362 505	Quick start manual
362 516	CD user manual Gymna devices multi language





Optional accessories

380 439	Carrying bag touchscreen electro devices
341 099	Contact gel, 5L
341 121	Pump, 5L

	COMBI 400V & 400 VIP	DUO 400V & 400 VIP	COMBI 400 & 400M	DUO 400 & 400M	PULSON 400 & 400 M
Therapies					
Electrotherapy (2 independent channels)	■	■	■	■	■
Ultrasound therapy (1 & 3 MHz)	■	■	■	■	■
Laser therapy (optional)	■	■	■	■	■
Combination therapy	■	■	■	■	■
Simultaneous therapy	■	■	■	■	■
Vacuum	■	■	■	■	■
User-interface					
Full colour TFT display, 10.4 inch (SVGA: 800 X 600 pixels)	■	■	■	■	■
Touch screen	■	■	■	■	■
Customisation wizard	■	■	■	■	■
Colour guided therapy methods	■	■	■	■	■
Enlarged therapy windows in dashboard design	■	■	■	■	■
2 separate intensity regulators	■	■	■	■	■
Guided Therapy System (GTS)	■	■	■	■	■
Medical E-book: anatomical library	■	■	■	■	■
Help and Clinical information screens	■	■	■	■	■
Direct therapy keys	■	■	■	■	■
Protocols: objectives, list of indications, selection for each body area	■	■	■	■	■
Protocols: cellular effects [heal the tissue]	■	■	■	■	■
3D pictures of electrode placement	■	■	■	■	■
Diagnostics (S-D curve, Rheobase, Chronaxy, ...)	■	■	■	■	■
Contra-indications list	■	■	■	■	■
Memory (free locations)					
500 for favourites/own programs	■	■	■	■	■
200 for diagnostic results	■	■	■	■	■
100 for own sequential programs	■	■	■	■	■
50 for shared programs on multiple devices	■	■	■	■	■

		COMBI 400V & 400 VIP	DUO 400V & 400 VIP	COMBI 400 & 400M	DUO 400 & 400M	PULSON 400 & 400 M
Electrotherapy 						
Unidirectional currents allowed with: 						
Direct current		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Rectangular pulse 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
2-5 current (Ultra Reiz) 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Triangular pulse 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
MF rectangular pulse 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Iontophoresis- MF rectangular pulse 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Iontophoresis- direct current 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Diadynamic currents						
MF 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
DF 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
RS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
CP 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
LP 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
TENS currents						
Conventional TENS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Low frequency TENS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Burst TENS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
High frequency TENS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Random Frequency TENS 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Han Stim (via painrelief)		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
NMES currents						
Rectangular surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Triangular surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Biphasic surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Intrapuls interval surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Russian stimulation 		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
2-pole MF surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Isoplanar vector field surge		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

		COMBI 400V & 400 VIP	DUO 400V & 400 VIP	COMBI 400 & 400M	DUO 400 & 400M	PULSON 400 & 400 M
Interferential currents						
2-pole Medium Frequency			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Isoplanar vector field			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Dipole vector field			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Classical interferential			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Micro current						
Micro current			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Micro current modulated			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Micro current surge			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
High Voltage (HVPC)						
High Voltage			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
High Voltage surge			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Diagnostic programs						
Rheobase and Chronaxy			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Rheobase and AQ			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
S-D curve rectangular			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
S-D curve triangular			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
S-D curve rectangular + triangular			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Pain points			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Diagnose stress fracture			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Iontophoresis programs			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Phonophoresis programs			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Constant Voltage / Constant Current			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Ultrasound therapy						
Hybrid treatment head 4 cm ² (1 & 3 MHz, multifrequent)			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Hybrid treatment head 1 cm ² (1 & 3 MHz, multifrequent)			<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

	COMBI 400V & 400 VIP	DUO 400V & 400 VIP	COMBI 400 & 400M	DUO 400 & 400M	PULSON 400 & 400 M
Laser therapy 					
Monoprobe 400: max. average power: 70,5 mW	○		○		
Clusterprobe 400: max. average power: 4 x 12,6 mW	○		○		
Combination therapy 					
See electro currents with 	■		■		
Simultaneous therapy					
Electrotherapy (2-pole & 4-pole) + Laser (optional)	■		■		
Electrotherapy (2-pole & 4-pole) + Ultrasound	■		■		
Ultrasound + Laser (optional)	■		■		
Electrotherapy (2-pole) + Electrotherapy (2-pole)	■	■	■	■	
Vacuum 					
2 independent channels	■	■			
Electronic vacuum control	■	■			
Continuous & pulsed rhythm	■	■			
Massage effect	■	■			
Control via Combi 400 or Duo 400	■	■			
Vacuum screen in dashboard design	■	■			

■ = Standard
○ = Optional