

Gymna ShockMaster 500

Treat under the best possible conditions

Product data sheet

User friendly

- Two output channels for simple and rapid switch between hand pieces



- Noiseless and vibration-free oil compressor, easy accessible
- Soft-close drawer with shock absorbers
- Unique storage for the applicators in the drawer



- Rapid software update via USB port



Always up to date


Shockwave therapy is evolving constantly. Thanks to treatment software updates, you stay abreast of new developments:

- **Unique! patient Guided Therapy System:** the intelligent Shockmaster patient Guided Therapy System will automatically adapt treatment level, based on your patient's condition
- **Total therapy plan:** every application is divided in different steps: pGTS, tendon treatment, trigger point treatment, muscle treatment with D-Actor, muscle treatment with V-Actor
- Treatment parameters are adapted in function of the latest information
- New applications are added, as well as relevant treatment protocols




Intuitive user interface

- **Easy application search:** submenu's based on the affected tissue

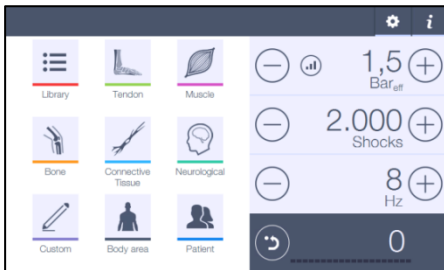

- **Unique soft start function** gradual pressure buildup for higher comfort for patient
- **Unique anatomical library:** function, origin, insertion, innervation and blood supply for all main muscles
- **Unique Marked treatment area and referred pain with trigger point**

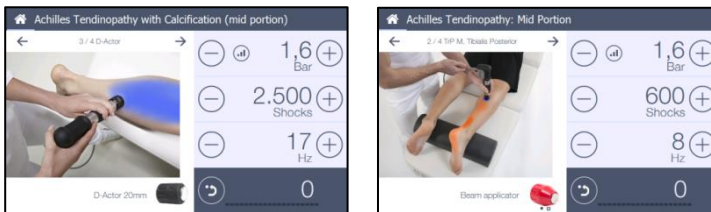

- Medical information on indications treated by shockwaves
- Patient follow up system via VAS (pain score)



- **Safety restriction:** preventing the use of an inappropriate handpiece
- **Pop-up** after 1.25 million shocks, when revision is needed

NEW Handpiece Short

- Ergonomic design: lightweight and well balanced
- Integrated soft grip materials
- Extended energy range: starting level at 0,3bar_{eff}
- Broader application: more patients, more applications
- Compatible with new Fascia and Spine applicators
- Backwards Compatibility: the Handpiece Short may be used on all ShockMaster devices provided that the mandatory software update is performed



Fascia applicator set



Spine applicator set

Gymna Club

Free subscription on the Gymna club: Extra information about shockwave therapy with the Gymna ShockMaster: literature, technical support, software updates for free, extra marketing material

Technical specifications

Characteristics

Dimensions (H x D x W cm)	102 x 58 x 57
Weight (kg)	65
Integrated trolley	4 extra-large swiveling wheels for superior mobility Made in aluminum Side protection for the device
Compressor	In trolley (oil), noiseless, vibration free 8 bar
Connections	2
Touch screen	7inch Full color 180° view

User interface

Plug and play	Choose your own settings
List of indications	43 indications
Total therapy plan	Every indication divided in different steps
Pre-programmed settings	> 90 programs
Body part	Choose indications via a body area
Marked treatment area	Treatment area marked on the pictures during treatment
Trigger point referred pain indication	
Quick start	Tendon applications Muscle applications Bone applications Connective applications Neurological applications Custom applications
Patient database	Save up to 2500 patients Visual analogue pain scale Follow up pain evolution
Backup system	On usb stick
Medical e-book	Anatomical library (muscles, nerves and blood vessels) Information about shockwave indications Information about applicators Contra-indications
Pre-set counter	Supported via illuminated bar
Soft start	Slowly building up intensity automatically
Software update	Via USB port

Performance

Frequency range	0,5 - 21 Hz (in steps of 0,1 Hz)
Max frequency with V-Actor	35 Hz
Pressure range	0,3 - 5 bar _{eff} (in steps of 0.1 bar _{eff})

Hand piece and applicators

On/off function	Electronic, no air cooling
Weight (kg)	0,48 kg
Ergonomics	Integrated soft rubber protection
V-Actor vibration therapy	Available
Revision kit (standard)	4 revisions

Applicators (standard)	Do it yourself revision Beam applicator D-Actor 20 mm applicator
Applicators (optional)	Deep Impact applicator 15mm applicator 6mm applicator D-Actor 35 mm applicator Focus lens applicator

Marketing material

Waiting room Poster
Certificate

Support

Gymna Club Scientific literature, Technical support, Free software updates,
marketing material, letter for the doctor

Ordering information

<u>Article</u>	<u>Article number</u>
Gymna ShockMaster 500	342.500

Standard included:

Device ShockMaster 500, Trolley ShockMaster 500, RSW Handpiece Short, Beam applicator, Deep Impact applicator, 20mm D-Actor applicator, 4 revision sets (2 projectiles and tubes inside handpiece box + 1 revision kit HP Short), 1 anti-skid protection mat, 1 bottle 500 ml contact gel, 1 tissue box, 1 waiting room poster, 1 certificate

Accessories:

RSW Handpiece Short	342.716
Upgrade pack SCM500	342.522
Handpiece short, Body replacement set	342.749
Revision Kit	309.771
Revision kit Handpiece Short	342.705
ShockMaster V-Actor II hand piece	338.239
Fascia applicator set (HP Short)	342.727
Spine applicator set (HP Short)	342.738
15mm Applicator	338.613
Beam Applicator	338.618
Deep Impact Applicator	338.624
Focus Lens	338.630
Acupuncture 6 mm	338.602
D-ACTOR 20mm S	338.646
D-ACTOR 35mm S	338.657
Contact 500 ml	341.088
Contact 5 L	341.099
Protective anti-skid mat	338.217
V-ACTOR 25mm	338.668
V-ACTOR 40mm	338.679
ShockMaster Poster	338.382
ShockMaster 500 Trolley	338.960