

Product datasheet

Tilting Advanced

Height-adjustable verticalisation table with two motors. One motor ensures effortless and smooth tilting; the other adjusts height and facilitates patient transfer. The Tilting is the ideal standing aid for neuromotor revalidation and adaptation exercises from lying to standing. Another area where it can be applied is general orthostatic examination, more specifically examination and treatment of vasovagal syncope.

Standard features

- Height adjustment and tilting by electric motors
- Both motors operated by hand switch
- Footrests can be set individually and can be folded down completely
- Large casters (lockable) guarantee the best possible mobile use
- Fitted with a tilt angle indicator
- Patient hoist compatible, 14 cm clearance below chassis
- Optional: battery power and separate leg fixation
- Gymna Stability Profile (GSP)
- Cushions with 40 mm padding
- Cushions in a choice of 11 colours

Ordering data

121 095 Tilting Advanced with hand switch

Optional

125 606 Separate leg fixation

125 607 back-up battery power supply

Accessories

112 214 pelvic/thoracic strap

112 215 Knee strap


112 218 Work top

115 768 Leg fastening strap, set of 2

125 655 Armrest, set of 2

xxx xxx Terry-cloth cover, stretch material
(available in 9 colours 'Style')

Specifications

Length x width:	± 202 x 67cm
Height adjustment:	± 62 – 110 cm
Weight:	± 134 kg
Max. lifting capacity:	150 kg
Electric motors:	8.000 N
Elapsed time high/low:	25 s. (35 s. on battery)
Elapsed time tilting:	35 s. (45 s. on battery)
Power supply:	230 VAC, 50 Hz, 0,8 A 120 VAC, 60 Hz, 0,8 A
Insulation class:	II
Moister resistance:	IP 21
Classification:	I
Conformity:	IEC/EN 60601-1 IEC/EN 60601-2 MDD 93/42/EEC
	
Guarantee:	3 years (upholstery 1 year)



Gymna Stability Profile

Tilting Advanced




Optional arm rests



Optional work top

Functionality/dimensions

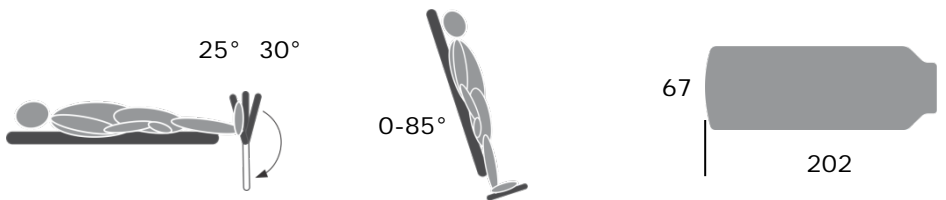

Cushion colours

Style



Fashion

